



Dual Bore Cannulas

MedOne coaxial dual bore cannulas allow injection of a fluid, such as perfluorocarbon or staining dye through a central lumen while simultaneously allowing pressure release via passive egress of fluid or air through a vent hole in the outer lumen. MedOne utilizes an extremely thin-wall inner cannula for an excellent flow rate, even for 25g cases. Dead space in the front of the cannula hub is filled to reduce air bubbles during injection. Inner cannula extends 3mm beyond outer cannula. Overall length 31mm.

- 3237 Dual Bore Cannula .9mm (20g)
- 3239 Dual Bore Cannula .6mm (23g)
- 3240 Dual Bore Cannula .5mm (25g)



MVR Knives

MedOne MVR knives are sharp and consistent, the two most important attributes of a surgical knife. The diamond architecture of the blade tip ensures a clean, precise incision every time. The 20g MVR creates a fluid tight seal for 20g instruments. The 23g MVR (3268/3269) can be used to create sclerotomies for use with 23g cannula systems but will *not fit through an entry cannula*.

- 3266 AdvantEdge® MVR Knife (20g)
- 3267 AdvantEdge® MVR Knife (20g) Angled
- 3268 AdvantEdge® MVR Knife (23g)
- 3269 AdvantEdge® MVR Knife (23g) Angled

- Available straight or angled for surgeon comfort
- Color coded handles for easy size identification
20g-Yellow handle, 23g-Blue handle



Narrow 23g MVR

Designed specifically to *fit through* a 23 gauge entry cannula and can be used as a knife or pick while using a 23 gauge cannula system.



3272 AdvantEdge® MVR Knife (23g) Narrow

Contact us for more information on our complete range of products for retina surgery

Exactly
what you're looking for




MedOne Surgical, Inc.
670 Tallevast Road, Sarasota, FL 34243 USA
TEL: 941.359.3129 Toll Free: 866.633.6631
Fax: 941.359.1708

www.MedOne.com



PolyTip® VFI Cannula 23g (10mm)

Similar to our popular 7mm VFI cannula, the 3241 has a 10mm cannula length to allow for improved access across the eye and visualization of the cannula tip during use. Designed for aspiration of 1000cs oil through a 23 gauge entry cannula.

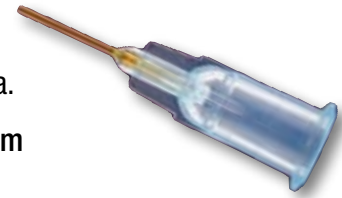

3241 PolyTip® VFI Cannula (23g)10mm

FlexTip Brush

Soft silicone brush tip allows gentle brushing and manipulation of the retina surface.



3215 FlexTip Brush .9mm (20g) 1mm
32mm cannula with 1mm brush tip

Backflush Instruments

We combined our most popular cannulas and our backflush handle into one convenient product. Put precise control of aspiration at the tip of your finger. Available with active or passive aspiration. FlexTip™ products have a 1mm silicone tip.

3280	Backflush Cannula	(20g)	Passive
3281	Backflush FlexTip™1mm	(20g)	Passive
3282	Backflush Brush	(20g)	Passive
3283	Backflush Cannula	(20g)	Active
3284	Backflush FlexTip™ 1mm	(20g)	Active
3285	Backflush Brush	(20g)	Active
3286	Backflush Cannula	(23g)	Passive
3287	Backflush FlexTip™ 1mm	(23g)	Passive
3288	Backflush Brush	(23g)	Passive
3289	Backflush Cannula	(23g)	Active
3290	Backflush FlexTip™ 1mm	(23g)	Active
3291	Backflush Brush	(23g)	Active
3292	Backflush Cannula	(25g)	Passive
3293	Backflush FlexTip™ 1mm	(25g)	Passive
3294	Backflush Brush	(25g)	Passive
3295	Backflush Cannula	(25g)	Active
3296	Backflush FlexTip™ 1mm	(25g)	Active
3297	Backflush Brush	(25g)	Active

